

Biosafety Cabinet Class II B2



Biosafety Cabinet Class II B2 is an equipment designed to provide a safe working environment for handling biological materials, such as microorganisms, viruses, and other potentially infectious agents

Airflow Pattern :

- Total Exhaust 100% exhaust system, where all air is exhausted through HEPA filters to the outside environment, minimizing the risk of recirculating potentially contaminated air.

Applications :


- **Handling Hazardous Agents :** These are substances that require a high level of containment, including certain pathogens and toxic chemicals.
- **Pharmaceutical and Biotechnology Research :** This involves research and development activities that necessitate strict containment and safety measures.
- **Clinical and Research Laboratories :** These facilities handle materials that pose significant risks to personnel and the environment.

Benefits :


- **Enhanced Safety :** Provides a high level of protection for both the operator and the environment by ensuring that all potentially contaminated air is exhausted.
- **Compliance :** Helps laboratories comply with stringent biosafety regulations and guidelines for handling hazardous materials.
- **Versatility :** Suitable for a wide range of applications requiring strict containment and safety measures.



Product Description			
Model	Nominal Width	Inner Chamber (WxDxH) mm	Outer Chamber (WxDxH) mm
AS-BSC-B2-4MS 4SS	4 Ft.	1120 x 640 x 640 mm	1260 x 800 x 2200 mm
AS-BSC-B2-6MS 6SS	6 Ft.	1740 x 650 x 640 mm	1860 x 800 x 2200 mm
Customize size on request			

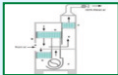

Total Exhaust System
No Air Recirculation

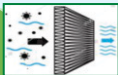

HEPA Filter
Dual HEPA Filtration for Inflow and Exhaust



Stainless Steel Interior
Easy to clean, Corrosion resistant



UV Germicidal Lamp
Optional for added sterilization



Microprocessor Control Panel
Intuitive and Programmable



Low Noise, High Efficiency
Quiet operation with powerful blower



Ettore Darvi Pvt. Ltd.
314, FIE Patparganj Industrial Area,
New Delhi-110092
E : sales@ettoredarvi.com
W : www.ettoredarvi.com